

**OneTaq® Quick-Load®  
2X Master Mix with  
Standard Buffer**



**M0486S**



**100 reactions (50 µl vol) Lot: 0121311**

**RECOMBINANT Store at -20°C Exp: 11/14**

**Description:** OneTaq Quick-Load 2X Master Mix with Standard Buffer is an optimized, ready-to-use blend of Taq and Deep Vent<sub>R</sub> DNA Polymerases ideally suited to routine PCR applications from templates, including pure DNA solutions, bacterial colonies and cDNA products. The 3'→5' exonuclease activity of Deep Vent DNA Polymerase increases the fidelity and robust amplification of Taq DNA Polymerase (1). The convenient quick-load master mix formulation contains dNTPs, MgCl<sub>2</sub>, buffer components and stabilizers as well as two commonly used tracking dyes for DNA gels. On a 1% agarose gel in 1X TBE, Xylene Cyanol FF migrates at ~4 kb and Tartrazine migrates at ~10 bp. Both dyes are present in concentrations that do not mask co-migrating DNA bands.

**Source:** An *E. coli* strain that carries the Taq DNA Polymerase gene from *Thermus aquaticus* YT-1 and an *E. coli* strain that carries the Deep Vent<sub>R</sub> DNA Polymerase gene from *Pyrococcus* species GB-D.

**Applications:**

- High Sensitivity PCR
- High Throughput PCR
- Routine PCR
- AT-rich PCR
- Colony PCR
- Long PCR (up to ~6 kb genomic)

**Reaction Conditions:** 1X OneTaq Quick-Load Master Mix with Standard Buffer, DNA template and primers in a total reaction volume of 50 µl.

**1X OneTaq Quick-Load Master Mix with Standard Buffer:**

- 20 mM Tris-HCl (pH 8.9 @ 25°C)
- 1.8 mM MgCl<sub>2</sub>
- 22 mM NH<sub>4</sub>Cl
- 22 mM KCl
- 0.2 mM dNTPs
- 5% glycerol
- 0.06% IGEPAL® CA-630
- 0.05% Tween® 20
- Xylene Cyanol FF
- Tartrazine
- 25 units/ml OneTaq DNA Polymerase

**Unit Definition:** One unit is defined as the amount of enzyme that will incorporate 15 nmol of dNTP into acid insoluble material in 30 minutes at 75°C.


**Unit Assay Conditions:** 1X ThermoPol® Reaction Buffer, 200 µM dNTPs including [<sup>3</sup>H]-dTTP and 200 µg/ml activated Calf Thymus DNA.

**Heat Inactivation:** No

**Quality Control Assays**

**5 kb Lambda PCR:** 25 cycles of PCR amplification of 5 ng Lambda DNA with 1X OneTaq Quick-Load Master Mix with Standard Buffer in a 25 µl reaction in the presence of 0.2 µM primers resulted in the expected 5 kb product.

**Note:** Product specifications for individual components in the OneTaq Quick-Load 2X Master Mix with Standard Buffer are available separately.

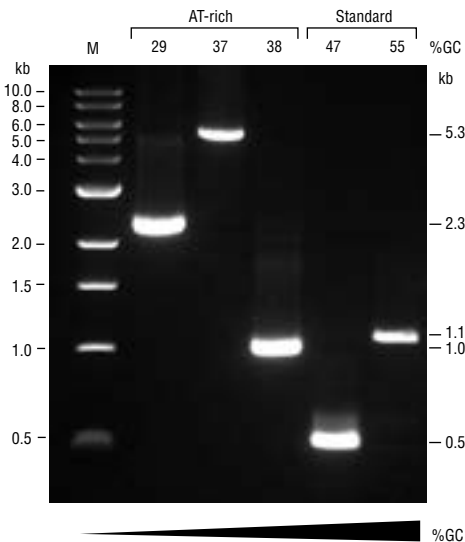

Amplification of a selection of sequences with varying GC content from human and *C. elegans* genomic DNA using OneTaq DNA Polymerase. Amplicon sizes are indicated next to gel, and GC content is indicated above gel. Marker M is the 1 kb DNA Ladder (NEB #N3232).

**PCR**

The Polymerase Chain Reaction (PCR) is a powerful and sensitive technique for DNA amplification (2). Taq DNA Polymerase is an enzyme widely used in PCR (3). The following guidelines are provided to ensure successful PCR using New England Biolabs' OneTaq Quick-Load 2X Master Mix with Standard Buffer. These guidelines cover routine PCR reactions. Amplification of templates with high GC content, high secondary structure or low template concentrations may require further optimization.

**Reaction Setup:**

We recommend assembling all reaction components on ice and quickly transferring the reactions to a thermocycler preheated to the denaturation temperature (94°C).

COMPONENT	25 µl REACTION	50 µl REACTION	FINAL CONCENTRATION
OneTaq Quick-Load 2X Master Mix with Standard Buffer	12.5 µl	25 µl	1X
10 µM Forward Primer	0.5 µl	1 µl	0.2 µM
10 µM Reverse Primer	0.5 µl	1 µl	0.2 µM
Template DNA	variable	variable	<1,000 ng
Nuclease-Free Water	to 25 µl	to 50 µl	

**Notes:** Gently mix the reaction. Collect all liquid to the bottom of the tube by a quick spin if necessary. Overlay the sample with mineral oil if using a PCR machine without a heated lid.

Transfer PCR tubes to a PCR machine and begin thermocycling:

**Thermocycling Conditions for a Routine PCR:**

STEP	TEMP	TIME
Initial Denaturation	94°C	30 seconds
30 Cycles	94°C	15–30 seconds
	45–68°C	15–60 seconds
Final Extension	68°C	1 minute/kb
Hold	68°C	5 minutes
	4–10°C	

**General Guidelines:**

1. **Template:**  
Use of high quality, purified DNA templates greatly enhances the success of PCR reactions. Recommended amounts of DNA template for a 50 µl reaction are as follows:

DNA	AMOUNT
Genomic	1 ng–1 µg
Plasmid or Viral	1 pg–1 ng

2. **Primers:**  
Oligonucleotide primers are generally 20–40 nucleotides in length and ideally have a GC content of 40–60%. Computer programs such as Primer3 (<http://frodo.wi.mit.edu/primer3>) can be used to design or analyze primers. The final concentration of each primer in a PCR reaction may be 0.05–1 µM, typically 0.2 µM.

3. **Mg<sup>++</sup> and Additives:**  
Mg<sup>++</sup> concentration of 1.5–2.0 mM is optimal for most PCR products generated with OneTaq DNA Polymerase. The final Mg<sup>++</sup> concentration in 1X OneTaq Quick-Load Master Mix with Standard Buffer is 1.8 mM. This supports satisfactory amplification of most amplicons. However, Mg<sup>++</sup> can be further optimized in 0.2 mM increments using MgCl<sub>2</sub> (sold separately).

For amplification of difficult targets, like GC-rich sequences, we recommend OneTaq 2X Master Mix with GC Buffer (sold separately). Alternatively, DMSO (4) or formamide (5) may be used.

4. **Denaturation:**  
An initial denaturation of 30 seconds at 94°C is sufficient to amplify most targets from pure DNA templates. For difficult templates such as GC-rich sequences, a longer denaturation of 2–4 minutes at 94°C is recommended prior to PCR cycling to fully denature the template. Alternatively, use OneTaq Quick-Load 2X Master Mix with GC Buffer. With colony PCR, an initial 2–5 minute denaturation at 94°C is recommended to lyse cells.

During thermocycling a 15–30 second denaturation at 94°C is recommended.

(see other side)

5. **Annealing:**  
The annealing step is typically 15–60 seconds. Annealing temperature is based on the  $T_m$  of the primer pair and is typically 45–68°C. Annealing temperatures can be optimized by doing a temperature gradient PCR starting 5°C below the calculated  $T_m$ . We recommend using NEB's  $T_m$  Calculator, available at [www.neb.com/TmCalculator](http://www.neb.com/TmCalculator) to determine appropriate annealing temperatures for PCR.

6. **Extension:**  
The recommended extension temperature is 68°C. Extension times are generally 1 minute per kb. A final extension of 5 minutes at 68°C is recommended.

7. **Cycle Number:**  
Generally, 25–35 cycles yields sufficient product. Up to 45 cycles may be required to detect low copy number targets.

8. **2-step PCR:**  
When primers with annealing temperatures of 68°C or above are used, a 2-step thermocycling protocol (combining annealing and extension into one step) is possible.

9. **PCR Product:**  
The majority of the PCR products generated using *OneTaq* DNA Polymerase contain dA overhangs at the 3' end; therefore the PCR products can be ligated to dT/dU-overhang vectors.

**Notes:**

*OneTaq* Quick-Load 2X Master Mix with Standard Buffer is stable for fifteen freeze-thaw cycles when stored at –20°C.

*OneTaq* Quick-Load 2X Master Mix with Standard Buffer is also stable for one month at 4°C, so for frequent use, an aliquot may be kept at 4°C.

**References:**

1. Barnes, W.M. (1994) *Proc. Natl. Acad. Sci. USA*, 91, 2216–2220.
2. Saiki R.K. et al. (1985) *Science*, 230, 1350–1354.
3. Powell, L.M. et al. (1987) *Cell*, 50, 831–840.
4. Sun, Y., Hegamyer, G. and Colburn, N. (1993) *Biotechniques*, 15, 372–374.
5. Sarkar, G., Kapelner, S. and Sommer, S.S. (1990) *Nucleic Acids Res.* 18, 7465.

**Companion Products Sold Separately:**

Magnesium Chloride (MgCl<sub>2</sub>) Solution  
#B9021S 6.0 ml

*OneTaq* Quick-Load 2X Master Mix with GC Buffer  
#M0487S 100 Reactions  
#M0487L 500 Reactions

*OneTaq* Hot Start Quick-Load 2X Master Mix with Standard Buffer  
#M0488S 100 Reactions  
#M0488L 500 Reactions

*OneTaq* Hot Start Quick-Load 2X Master Mix with GC Buffer  
#M0489S 100 Reactions  
#M0489L 500 Reactions

*OneTaq* DNA Polymerase  
#M0480S 200 units  
#M0480L 1,000 units  
#M0480X 5,000 units

*OneTaq* Hot Start DNA Polymerase  
#M0481S 200 units  
#M0481L 1,000 units  
#M0481X 5,000 units



Limited Label License: This product is covered by certain patents outside the United States (the "Non-US Patents") corresponding to the expired U.S. Patent No. 5,436,149 (claims 6-16). The purchase of this product conveys to the purchaser only the limited, non-transferable right under the Non-US Patents to use the purchased quantity of the product for the purchaser's own research by the purchaser only. No rights are granted to the purchaser with respect to the Non-US Patents to sell, modify for resale or otherwise transfer this product, either alone or as a component of another product, to any third party. Takara Bio reserves all other rights to the Non-US Patents, and this product may not be used in any manner other than as provided herein. For information on obtaining a license with respect to the Non-US Patents to use this product for purposes other than research, please contact Takara Bio Inc., Seta 3-4-1, Otsu, Shiga 520-2193, Japan (Fax +81-77-543-9254).

NEW ENGLAND BIOLABS®, ONETAQ®, QUICK-LOAD® and THERMOPOL® are registered trademarks of New England Biolabs, Inc. DEEP VENT™ is a trademark of New England Biolabs, Inc.

IGEPAL® is a registered trademark of Rhodia Operations. TWEEN® is a registered trademark of Uniqema Americas LLC.