

β 1-3 Galactosidase



1-800-632-7799
info@neb.com
www.neb.com



P0726S 007150817081

P0726S



500 units 10,000 U/ml Lot: 0071508

RECOMBINANT Store at -20°C Exp: 8/17

Description: β 1-3 Galactosidase is a highly specific exoglycosidase that catalyzes the hydrolysis of β 1-3 and, at a much lower rate, β 1-6 linked D-galactopyranosyl residues from oligosaccharides. The approximate kinetic data show > 100-fold preference for β 1-3 over β 1-6 linkages (1,2) and > 500-fold preference from β 1-3 over β 1-4 linkages (3).

New Reaction Buffer

β 1-3 Galactosidase



1-800-632-7799
info@neb.com
www.neb.com



P0726S 007150817081

P0726S



500 units 10,000 U/ml Lot: 0071508

RECOMBINANT Store at -20°C Exp: 8/17

Description: β 1-3 Galactosidase is a highly specific exoglycosidase that catalyzes the hydrolysis of β 1-3 and, at a much lower rate, β 1-6 linked D-galactopyranosyl residues from oligosaccharides. The approximate kinetic data show > 100-fold preference for β 1-3 over β 1-6 linkages (1,2) and > 500-fold preference from β 1-3 over β 1-4 linkages (3).


New Reaction Buffer

Specificity:

↓
Gal β 1 – 3 R
> β 1 – 6 R
>> β 1 – 4 R

Detailed Specificity: The GlcNAc(β 1-6) residue is the only anomeric configuration that can effect the specificity of the enzyme enabling cleavage of the non-reducing β 1-4Galactose (Fig. 1).

A) 0.1 nm/ μ l substrate, 10 units, 1 hr incubation

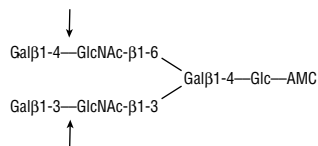

B) 0.1 nm/ μ l substrate, 0.625 units, 1 hr incubation

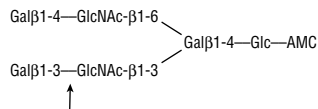



Figure 1: Selling concentration of the enzyme will cut the β 1-4Galactose linkage as shown in (A) due to the adjacent GlcNAc β 1-6 anomer. This cleavage will not occur if the selling concentration of the enzyme is diluted 16-fold, as shown in (B).

Specificity:

↓
Gal β 1 – 3 R
> β 1 – 6 R
>> β 1 – 4 R

Detailed Specificity: The GlcNAc(β 1-6) residue is the only anomeric configuration that can effect the specificity of the enzyme enabling cleavage of the non-reducing β 1-4Galactose (Fig. 1).

A) 0.1 nm/ μ l substrate, 10 units, 1 hr incubation

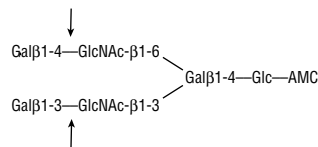

B) 0.1 nm/ μ l substrate, 0.625 units, 1 hr incubation

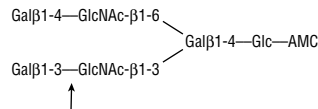


Figure 1: Selling concentration of the enzyme will cut the β 1-4Galactose linkage as shown in (A) due to the adjacent GlcNAc β 1-6 anomer. This cleavage will not occur if the selling concentration of the enzyme is diluted 16-fold, as shown in (B).

Source: Cloned from *Xanthomonas manihotis* and expressed in *E. coli* (4).

Supplied in: 50 mM NaCl, 20 mM Tris-HCl (pH 7.5 @ 25°C) and 1 mM EDTA.

Reagents Supplied with Enzyme:

10X GlycoBuffer 1
100X BSA

Reaction Conditions:

1X GlycoBuffer 1:
50 mM Sodium Acetate (pH 5.5 @ 25°C)
and 5 mM CaCl_2 . Supplement with 100 $\mu\text{g}/\text{ml}$ BSA.
Incubate at 37°C .

Optimal incubation times and enzyme concentrations must be determined empirically for a particular substrate.

Unit Definition: One unit is defined as the amount of enzyme required to cleave > 95% of the terminal β -D-galactose from 1 nmol of Gal β 1-3GlcNAc β 1-3Gal β 1-4Glc-7-amino-4-methyl-coumarin (AMC), in 1 hour at 37°C in a total reaction volume of 10 μl .

Source: Cloned from *Xanthomonas manihotis* and expressed in *E. coli* (4).

Supplied in: 50 mM NaCl, 20 mM Tris-HCl (pH 7.5 @ 25°C) and 1 mM EDTA.

Reagents Supplied with Enzyme:

10X GlycoBuffer 1
100X BSA

Reaction Conditions:

1X GlycoBuffer 1:
50 mM Sodium Acetate (pH 5.5 @ 25°C)
and 5 mM CaCl_2 . Supplement with 100 $\mu\text{g}/\text{ml}$ BSA.
Incubate at 37°C .

Optimal incubation times and enzyme concentrations must be determined empirically for a particular substrate.

Unit Definition: One unit is defined as the amount of enzyme required to cleave > 95% of the terminal β -D-galactose from 1 nmol of Gal β 1-3GlcNAc β 1-3Gal β 1-4Glc-7-amino-4-methyl-coumarin (AMC), in 1 hour at 37°C in a total reaction volume of 10 μl .

Unit Definition Assay: Two fold dilutions of β 1-3 Galactosidase are incubated with 1 nmol AMC-labeled substrate in 1X GlycoBuffer 1, supplemented with 100 $\mu\text{g}/\text{ml}$ BSA, in a 10 μl reaction. The reaction mix is incubated for 1 hour at 37°C . Separation of reaction products are visualized via thin layer chromatography (1).

Specific Activity: 17,000 units/mg

Molecular Weight: 66,000 daltons.

Quality Assurance: No contaminating exoglycosidase or proteolytic activity could be detected.

Quality Control Assays

Glycosidase Assay: 100 units of β 1-3 Galactosidase were incubated with 0.1 mM of fluorescently-labeled oligosaccharides and glycopeptides, in a 10 μl reaction for 20 hours at 37°C . The reaction products were analyzed by TLC for digestion of substrate.

(see other side)

CERTIFICATE OF ANALYSIS

Unit Definition Assay: Two fold dilutions of β 1-3 Galactosidase are incubated with 1 nmol AMC-labeled substrate in 1X GlycoBuffer 1, supplemented with 100 $\mu\text{g}/\text{ml}$ BSA, in a 10 μl reaction. The reaction mix is incubated for 1 hour at 37°C . Separation of reaction products are visualized via thin layer chromatography (1).

Specific Activity: 17,000 units/mg

Molecular Weight: 66,000 daltons.

Quality Assurance: No contaminating exoglycosidase or proteolytic activity could be detected.

Quality Control Assays

Glycosidase Assay: 100 units of β 1-3 Galactosidase were incubated with 0.1 mM of fluorescently-labeled oligosaccharides and glycopeptides, in a 10 μl reaction for 20 hours at 37°C . The reaction products were analyzed by TLC for digestion of substrate.

(see other side)

CERTIFICATE OF ANALYSIS

Physical Purity: Purified to > 95% homogeneity as determined by SDS-PAGE analysis using Coomassie Blue detection.

No other glycosidase activities were detected (ND) with the following substrates:

β -N-Acetyl-glucosaminidase:
GlcNAc β 1-4GlcNAc β 1-4GlcNAc-AMC ND

α -Fucosidase:
Fuc α 1-2Gal β 1-4Glc-AMCGal β 1-4
(Fuc α 1-3)GlcNAc β 1-3Gal β 1-4Glc-AMC ND

α -Galactosidase:
Gal α 1-3Gal β 1-4Gal-AMC ND

α -Neuraminidase:
Neu5Ac α 2-3Gal β 1-3GlcNAc β 1-3Gal β
1-4Glc-AMC ND

α -Mannosidase:
Man α 1-3Man β 1-4GlcNAc-AMC
Man α 1-6Man α 1-6(Man α 1-3)Man-AMC ND

Page 2 (P0726)

Physical Purity: Purified to > 95% homogeneity as determined by SDS-PAGE analysis using Coomassie Blue detection.

No other glycosidase activities were detected (ND) with the following substrates:

β -N-Acetyl-glucosaminidase:
GlcNAc β 1-4GlcNAc β 1-4GlcNAc-AMC ND

α -Fucosidase:
Fuc α 1-2Gal β 1-4Glc-AMCGal β 1-4
(Fuc α 1-3)GlcNAc β 1-3Gal β 1-4Glc-AMC ND

α -Galactosidase:
Gal α 1-3Gal β 1-4Gal-AMC ND

α -Neuraminidase:
Neu5Ac α 2-3Gal β 1-3GlcNAc β 1-3Gal β
1-4Glc-AMC ND

α -Mannosidase:
Man α 1-3Man β 1-4GlcNAc-AMC
Man α 1-6Man α 1-6(Man α 1-3)Man-AMC ND

Page 2 (P0726)

β -Glucosidase:
Glc β 1-4Glc β 1-4Glc-AMC ND

β -Xylosidase:
Xyl β 1-4Xyl β 1-4Xyl β 1-4Xyl-AMC ND

β -Mannosidase:
Man β 1-4Man β 1-4Man-AMC ND

Endo F₁, F₂, H:
Dansylated invertase high mannose. ND

Endo F₂, F₃:
Dansylated fibrinogen biantennary. ND

PNGase F:
Fluoresceinated fetuin triantennary. ND

Protease Assay: After incubation of 100 units of β 1-3 Galactosidase with 0.2 nmol of a standard mixture of proteins in a 20 μ l reaction, for 20 hours at 37°C, no proteolytic activity could be detected by SDS-PAGE.

β -Glucosidase:
Glc β 1-4Glc β 1-4Glc-AMC ND

β -Xylosidase:
Xyl β 1-4Xyl β 1-4Xyl β 1-4Xyl-AMC ND

β -Mannosidase:
Man β 1-4Man β 1-4Man-AMC ND

Endo F₁, F₂, H:
Dansylated invertase high mannose. ND

Endo F₂, F₃:
Dansylated fibrinogen biantennary. ND

PNGase F:
Fluoresceinated fetuin triantennary. ND

Protease Assay: After incubation of 100 units of β 1-3 Galactosidase with 0.2 nmol of a standard mixture of proteins in a 20 μ l reaction, for 20 hours at 37°C, no proteolytic activity could be detected by SDS-PAGE.

Note: Recommended storage temperature has changed to -20°C.

Avoid repeated freeze/thaw cycles.

References:

1. Wong-Madden, S.T. and Landry, D. (1995) *Glycobiology* 5, 19–28.
2. Guthrie, E.P. and Taron C. New England Biolabs, Inc., unpublished results.
3. Monks, B., New England Biolabs, Inc., unpublished results.
4. Taron, C.H. et al. (1995) *Glycobiology* 5, 603–610.



NEW ENGLAND BIOLABS® is a registered trademark of New England Biolabs, Inc.

This product is intended for research purposes only. This product is not intended to be used for therapeutic or diagnostic purposes in humans or animals.

Note: Recommended storage temperature has changed to -20°C.

Avoid repeated freeze/thaw cycles.

References:

1. Wong-Madden, S.T. and Landry, D. (1995) *Glycobiology* 5, 19–28.
2. Guthrie, E.P. and Taron C. New England Biolabs, Inc., unpublished results.
3. Monks, B., New England Biolabs, Inc., unpublished results.
4. Taron, C.H. et al. (1995) *Glycobiology* 5, 603–610.



NEW ENGLAND BIOLABS® is a registered trademark of New England Biolabs, Inc.

This product is intended for research purposes only. This product is not intended to be used for therapeutic or diagnostic purposes in humans or animals.