

# $\alpha$ 1-6 Mannosidase



1-800-632-7799  
info@neb.com  
www.neb.com



P0727S 009151216121

## P0727S



800 units    40,000 U/ml    Lot: 0091512  
RECOMBINANT    Store at 4°C    Exp: 12/16

**Description:**  $\alpha$ 1-6 Mannosidase is a highly specific exoglycosidase that removes unbranched  $\alpha$ 1-6 linked D-mannopyranosyl residues from oligosaccharides (1,2). When used in conjunction with  $\alpha$ 1-2,3 Mannosidase, the  $\alpha$ 1-6 Mannosidase will cleave  $\alpha$ 1-6 Mannose residues from branched carbohydrate substrates.

New Reaction Buffer

# $\alpha$ 1-6 Mannosidase



1-800-632-7799  
info@neb.com  
www.neb.com



P0727S 009151216121

## P0727S



800 units    40,000 U/ml    Lot: 0091512  
RECOMBINANT    Store at 4°C    Exp: 12/16

**Description:**  $\alpha$ 1-6 Mannosidase is a highly specific exoglycosidase that removes unbranched  $\alpha$ 1-6 linked D-mannopyranosyl residues from oligosaccharides (1,2). When used in conjunction with  $\alpha$ 1-2,3 Mannosidase, the  $\alpha$ 1-6 Mannosidase will cleave  $\alpha$ 1-6 Mannose residues from branched carbohydrate substrates.


New Reaction Buffer

### Specificity:


Man  $\alpha$  1 - 6 R

**Detailed Specificity:** Specificity can vary depending on incubation time and concentration of substrate (Figure 1).


#### A. 0.1 nm/ $\mu$ l substrate, 4 hour incubation

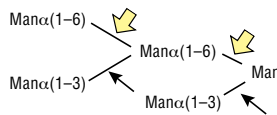

#### B. 0.1 nm/ $\mu$ l substrate, 4 hour incubation

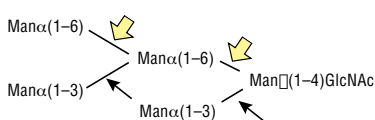

#### C. 0.1 nm/ $\mu$ l substrate, 18 hour incubation

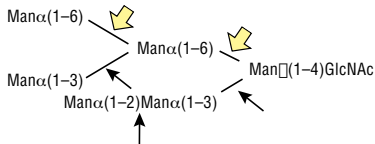


### Specificity:


Man  $\alpha$  1 - 6 R

**Detailed Specificity:** Specificity can vary depending on incubation time and concentration of substrate (Figure 1).


#### A. 0.1 nm/ $\mu$ l substrate, 4 hour incubation

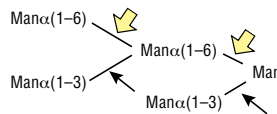


#### B. 0.1 nm/ $\mu$ l substrate, 4 hour incubation

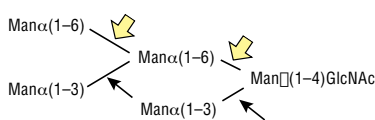


#### C. 0.1 nm/ $\mu$ l substrate, 18 hour incubation

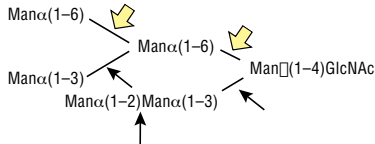

#### D. 0.05 nm/ $\mu$ l substrate, 18 hour incubation

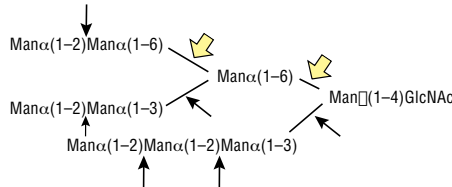


#### E. 0.045 nm/ $\mu$ l substrate, 18 hour incubation

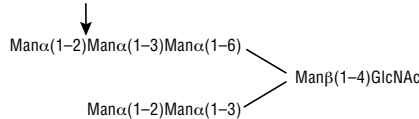


**Figure 1: Detailed specificity of  $\alpha$ 1-6 Mannosidase.** All reactions contained 32 units of  $\alpha$ 1-2,3 Mannosidase (NEB #P0729), 40 units of  $\alpha$ 1-6 Mannosidase, 1X GlycoBuffer 1 and 1X BSA in a total reaction volume of 10  $\mu$ l. Reactions were incubated at 37°C. The substrate depicted in (E) will not cut to completion. If this structure exists in any substrate it will be impervious to cleavage by  $\alpha$ 1-6 Mannosidase. Note: When used alone,  $\alpha$ 1-6 Mannosidase will still act only on linear substrates. When used in conjunction with  $\alpha$ 1-2,3 Mannosidase, the  $\alpha$ 1-6 Mannosidase will cleave  $\alpha$ 1-6 Mannose residues from branched carbohydrate substrates.

#### D. 0.05 nm/ $\mu$ l substrate, 18 hour incubation

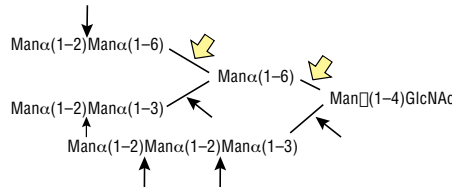

#### E. 0.045 nm/ $\mu$ l substrate, 18 hour incubation

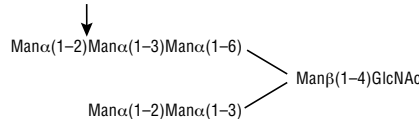

**Figure 1: Detailed specificity of  $\alpha$ 1-6 Mannosidase.** All reactions contained 32 units of  $\alpha$ 1-2,3 Mannosidase (NEB #P0729), 40 units of  $\alpha$ 1-6 Mannosidase, 1X GlycoBuffer 1 and 1X BSA in a total reaction volume of 10  $\mu$ l. Reactions were incubated at 37°C. The substrate depicted in (E) will not cut to completion. If this structure exists in any substrate it will be impervious to cleavage by  $\alpha$ 1-6 Mannosidase. Note: When used alone,  $\alpha$ 1-6 Mannosidase will still act only on linear substrates. When used in conjunction with  $\alpha$ 1-2,3 Mannosidase, the  $\alpha$ 1-6 Mannosidase will cleave  $\alpha$ 1-6 Mannose residues from branched carbohydrate substrates.

Note: *p*-nitrophenyl- $\alpha$ -D-mannopyranoside is NOT a substrate for this enzyme.

**Source:** Cloned from *Xanthomonas manihotis* and expressed in *E. coli* (2).

Supplied in: 50 mM NaCl, 20 mM Tris-HCl (pH 7.5 @ 25°C), 1 mM EDTA.

#### Reagents Supplied with Enzyme:

10X GlycoBuffer 1  
100X BSA

#### Reaction Conditions:

1X GlycoBuffer 1:  
50 mM Sodium Acetate (pH 5.5 @ 25°C)  
and 5 mM CaCl<sub>2</sub>. Supplement with 100  $\mu$ g/ml BSA. Incubate at 37°C.

**Unit Definition:** One unit is defined as the amount of enzyme required to cleave > 95% of the terminal  $\alpha$ -D-mannose from 1 nmol of Man $\alpha$ 1-6Man $\alpha$ 1-6Man-7-amino-4-methyl-coumarin (AMC), in 1 hour at 37°C in a total reaction volume of 10  $\mu$ l.

(see other side)

CERTIFICATE OF ANALYSIS

Note: *p*-nitrophenyl- $\alpha$ -D-mannopyranoside is NOT a substrate for this enzyme.

**Source:** Cloned from *Xanthomonas manihotis* and expressed in *E. coli* (2).

Supplied in: 50 mM NaCl, 20 mM Tris-HCl (pH 7.5 @ 25°C), 1 mM EDTA.

#### Reagents Supplied with Enzyme:

10X GlycoBuffer 1  
100X BSA

#### Reaction Conditions:

1X GlycoBuffer 1:  
50 mM Sodium Acetate (pH 5.5 @ 25°C)  
and 5 mM CaCl<sub>2</sub>. Supplement with 100  $\mu$ g/ml BSA. Incubate at 37°C.

**Unit Definition:** One unit is defined as the amount of enzyme required to cleave > 95% of the terminal  $\alpha$ -D-mannose from 1 nmol of Man $\alpha$ 1-6Man $\alpha$ 1-6Man-7-amino-4-methyl-coumarin (AMC), in 1 hour at 37°C in a total reaction volume of 10  $\mu$ l.

(see other side)

CERTIFICATE OF ANALYSIS

**Specific Activity:** ~ 137,000 units/mg

**Molecular Weight:** 51,000 daltons

**Unit Definition Assay:** Two fold dilutions of  $\alpha$ 1-6 Mannosidase are incubated with 1 nmol AMC-labeled substrate in 1X GlycoBuffer 1, supplemented with 100  $\mu$ g/ml BSA, in a 10  $\mu$ l reaction. The reaction mix is incubated for 1 hour at 37°C. Separation of reaction products are visualized via thin layer chromatography (1).

**Quality Assurance:** No contaminating exoglycosidase or proteolytic activity could be detected (ND).

### Quality Controls

**Glycosidase Assays:** 80 units of  $\alpha$ 1-6 Mannosidase were incubated with 0.1 mM of fluorescently-labeled oligosaccharides and glycopeptides, in a 10  $\mu$ l reaction for 20 hours at 37°C. The reaction products were analyzed by TLC for digestion of substrate.

Page 2 (P0727)

**Specific Activity:** ~ 137,000 units/mg

**Molecular Weight:** 51,000 daltons

**Unit Definition Assay:** Two fold dilutions of  $\alpha$ 1-6 Mannosidase are incubated with 1 nmol AMC-labeled substrate in 1X GlycoBuffer 1, supplemented with 100  $\mu$ g/ml BSA, in a 10  $\mu$ l reaction. The reaction mix is incubated for 1 hour at 37°C. Separation of reaction products are visualized via thin layer chromatography (1).

**Quality Assurance:** No contaminating exoglycosidase or proteolytic activity could be detected (ND).

### Quality Controls

**Glycosidase Assays:** 80 units of  $\alpha$ 1-6 Mannosidase were incubated with 0.1 mM of fluorescently-labeled oligosaccharides and glycopeptides, in a 10  $\mu$ l reaction for 20 hours at 37°C. The reaction products were analyzed by TLC for digestion of substrate.

Page 2 (P0727)

**$\beta$ -N-Acetylglucosaminidase:**  
GlcNAc $\beta$ 1-4GlcNAc $\beta$ 1-4GlcNAc-AMC ND

**$\alpha$ -N-Acetylgalactosaminidase:**  
GalNAc $\alpha$ 1-3(Fuc $\alpha$ 1-2)Gal $\beta$ 1-4Glc-AMC ND

**$\alpha$ -Fucosidase:**  
Gal $\beta$ 1-4(Fuc $\alpha$ 1-3)GlcNAc $\beta$ 1-3Gal $\beta$ 1-4Glc-AMC ND  
Fuc $\alpha$ 1-2Gal $\beta$ 1-4Glc-AMC ND

**$\alpha$ -Galactosidase:**  
Gal $\alpha$ 1-3Gal $\beta$ 1-4Gal-AMC ND

**$\beta$ -Galactosidase:**  
Gal $\beta$ 1-3GlcNAc $\beta$ 1-4Gal $\beta$ 1-4Glc-AMC ND  
Gal $\beta$ 1-4GlcNAc $\beta$ 1-3Gal $\beta$ 1-4Glc-AMC ND

**$\alpha$ -Neuraminidase:**  
Neu5Ac $\alpha$ 2-3Gal $\beta$ 1-3GlcNAc $\beta$ 1-3Gal $\beta$ 1-4Glc-AMC ND

**$\alpha$ -Mannosidase:**  
Man $\alpha$ 1-3Man $\beta$ 1-4GlcNAc-AMC ND

**$\beta$ -N-Acetylglucosaminidase:**  
GlcNAc $\beta$ 1-4GlcNAc $\beta$ 1-4GlcNAc-AMC ND

**$\alpha$ -N-Acetylgalactosaminidase:**  
GalNAc $\alpha$ 1-3(Fuc $\alpha$ 1-2)Gal $\beta$ 1-4Glc-AMC ND

**$\alpha$ -Fucosidase:**  
Gal $\beta$ 1-4(Fuc $\alpha$ 1-3)GlcNAc $\beta$ 1-3Gal $\beta$ 1-4Glc-AMC ND  
Fuc $\alpha$ 1-2Gal $\beta$ 1-4Glc-AMC ND

**$\alpha$ -Galactosidase:**  
Gal $\alpha$ 1-3Gal $\beta$ 1-4Gal-AMC ND

**$\beta$ -Galactosidase:**  
Gal $\beta$ 1-3GlcNAc $\beta$ 1-4Gal $\beta$ 1-4Glc-AMC ND  
Gal $\beta$ 1-4GlcNAc $\beta$ 1-3Gal $\beta$ 1-4Glc-AMC ND

**$\alpha$ -Neuraminidase:**  
Neu5Ac $\alpha$ 2-3Gal $\beta$ 1-3GlcNAc $\beta$ 1-3Gal $\beta$ 1-4Glc-AMC ND

**$\alpha$ -Mannosidase:**  
Man $\alpha$ 1-3Man $\beta$ 1-4GlcNAc-AMC ND

**$\beta$ -Glucosidase:**  
Glc $\beta$ 1-4Glc $\beta$ 1-4Glc-AMC ND

**$\alpha$ -Glucosidase:**  
Glc $\alpha$ 1-6Glc $\alpha$ 1-4Glc-AMC ND

**$\beta$ -Xylosidase:**  
Xyl $\beta$ 1-4Xyl $\beta$ 1-4Xyl $\beta$ 1-4Xyl-AMC ND

**$\beta$ -Mannosidase:**  
Man $\beta$ 1-4Man $\beta$ 1-4Man-AMC ND

**Endo F<sub>1</sub>, F<sub>2</sub>, H:**  
Dansylated invertase high mannose. ND

**Endo F<sub>2</sub>, F<sub>3</sub>:**  
Dansylated fibrinogen biantennary. ND

**PNase F:**  
Fluoresceinated fetuin triantennary. ND

**Protease Assay:** After incubation of 560 units of  $\alpha$ 1-6 Mannosidase with 0.2 nmol of a standard mixture of proteins for 20 hours at 37°C, no proteolytic activity could be detected by SDS-PAGE.

**$\beta$ -Glucosidase:**  
Glc $\beta$ 1-4Glc $\beta$ 1-4Glc-AMC ND

**$\alpha$ -Glucosidase:**  
Glc $\alpha$ 1-6Glc $\alpha$ 1-4Glc-AMC ND

**$\beta$ -Xylosidase:**  
Xyl $\beta$ 1-4Xyl $\beta$ 1-4Xyl $\beta$ 1-4Xyl-AMC ND

**$\beta$ -Mannosidase:**  
Man $\beta$ 1-4Man $\beta$ 1-4Man-AMC ND

**Endo F<sub>1</sub>, F<sub>2</sub>, H:**  
Dansylated invertase high mannose. ND

**Endo F<sub>2</sub>, F<sub>3</sub>:**  
Dansylated fibrinogen biantennary. ND

**PNase F:**  
Fluoresceinated fetuin triantennary. ND

**Protease Assay:** After incubation of 560 units of  $\alpha$ 1-6 Mannosidase with 0.2 nmol of a standard mixture of proteins for 20 hours at 37°C, no proteolytic activity could be detected by SDS-PAGE.

### References:

1. Wong-Madden, S.T. and Landry, D. (1995) *Glycobiology* 5, 19–28.
2. Guthrie, E. P., Taron, C.H., New England Biolabs, Inc., unpublished results.



NEW ENGLAND BIOLABS® is a registered trademark of New England Biolabs, Inc.

This product is intended for research purposes only. This product is not intended to be used for therapeutic or diagnostic purposes in humans or animals.

### References:

1. Wong-Madden, S.T. and Landry, D. (1995) *Glycobiology* 5, 19–28.
2. Guthrie, E. P., Taron, C.H., New England Biolabs, Inc., unpublished results.



NEW ENGLAND BIOLABS® is a registered trademark of New England Biolabs, Inc.

This product is intended for research purposes only. This product is not intended to be used for therapeutic or diagnostic purposes in humans or animals.