240 County Road Ipswich, MA 01938-2723

Tel 978-927-5054 Fax 978-921-1350 www.neb.com info@neb.com

New England Biolabs Product Specification

Product Name: Vent® (exo-) DNA Polymerase

Catalog #: M0257S/L
Concentration: 2,000 units/ml

Unit Definition: One unit is defined as the amount of enzyme that will incorporate 10 nmol of dNTP into acid-insoluble material in 30

minutes at 75°C

Shelf Life: 24 months
Storage Temp: -20°C

Storage Conditions: 10 mM Tris-HCl, 100 mM KCl, 1 mM DTT, 0.1 mM EDTA, 0.1 % Triton®X-100, 50 % Glycerol, (pH 7.4 @

25°C)

Specification Version: PS-M0257S/L v1.0

Effective Date: 02 Dec 2015

Assay Name/Specification (minimum release criteria)

Endonuclease Activity (Nicking) - A 50 μ l reaction in ThermoPol® Reaction Buffer containing 1 μ g of supercoiled PhiX174 DNA and a minimum of 200 units of Vent® (exo-) DNA Polymerase incubated for 4 hours at either 37°C or 75°C results in <10% conversion to the nicked form as determined by agarose gel electrophoresis.

Exonuclease Activity (Radioactivity Release) - A 50 μ l reaction in ThermoPol® Reaction Buffer containing 1 μ g of a mixture of single and double-stranded [3 H] *E. coli* DNA and a minimum of 200 units of Vent® (exo-) DNA Polymerase incubated for 4 hours at either 37° C or 75°C releases <0.1% of the total radioactivity.

Non-Specific DNase Activity (16 Hour) - A 50 μ l reaction in NEBuffer 2 containing 1 μ g of T3 DNA in addition to a reaction containing Lambda-HindIII DNA and a minimum of 2 units of Vent® (exo-) DNA Polymerase incubated for 16 hours at 37°C results in a DNA pattern free of detectable nuclease degradation as determined by agarose gel electrophoresis.

PCR Amplification (2.0 kb Lambda DNA) - A 25 μ l reaction in ThermoPol® Reaction Buffer in the presence of 200 μ M dNTPs and 0.2 μ M primers containing 5 ng Lambda DNA with 0.5 units of Vent® (exo-) DNA Polymerase for 30 cycles of PCR amplification results in the expected 2.0 kb product.


Phosphatase Activity (pNPP) - A 200 μ l reaction in 1M Diethanolamine, pH 9.8, 0.5 mM MgCl₂ containing 2.5 mM p-Nitrophenyl Phosphate (pNPP) and a minimum of 100 units Vent® (exo-) DNA Polymerase incubated for 4 hours at 37°C yields <0.0001 unit of alkaline phosphatase activity as determined by spectrophotometric analysis.

Protein Purity Assay (SDS-PAGE) - Vent® (exo-) DNA Polymerase is ≥ 95% pure as determined by SDS-PAGE analysis using Coomassie Blue detection.







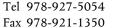

www.neb.com info@neb.com



New England Biolabs Product Specification

Assay Name/Specification (minimum release criteria)

qPCR DNA Contamination (*E. coli* Genomic) - A minimum of 2 units of Vent® (exo-) DNA Polymerase is screened for the presence of *E. coli* genomic DNA using SYBR® Green qPCR with primers specific for the *E. coli* 16S rRNA locus. Results are quantified using a standard curve generated from purified *E. coli* genomic DNA. The measured level of *E. coli* genomic DNA contamination is ≤ 1 *E. coli* genome.

RNase Activity (Extended Digestion) - A 10 μ l reaction in NEBuffer 4 containing 40 ng of a 300 base single-stranded RNA and a minimum of 1 μ l of Vent® (exo-) DNA Polymerase is incubated at 37°C. After incubation for 16 hours, >90% of the substrate RNA remains intact as determined by gel electrophoresis using fluorescent detection.

Single Stranded DNase Activity (FAM-Labeled Oligo) - A 20 μ l reaction in ThermoPol® Reaction Buffer containing a 10 nM solution of a fluorescent internal labeled oligonucleotide and a minimum of 20 units of Vent® (exo-) DNA Polymerase incubated for 30 minutes at either 37°C or 75°C yields <10% degradation as determined by capillary electrophoresis.

Derek Robinson

Director of Quality Control







Date

02 Dec 2015